

# OKAY<sup>®</sup> II & GM-30

## Y-Connector & Inflation Device

### YOUR LINK BETWEEN PATIENT AND CARE

Nipro Medical Corporation is a subsidiary of Nipro Corporation Japan, established in 1954 and a world-leader in healthcare with over 30,000 employees around the globe. We are committed to delivering innovative, high quality solutions in the form of renal, medical-surgical, and vascular products. This fundamental objective drives us to develop cutting-edge, value-added products for healthcare professionals that are safe, effective, and easy to use.

As an integral part of a vast global manufacturing network, Nipro Medical Corporation's North America division is responsible for managing the sales, marketing, and business operations of the United States and Canada – unified in our commitment to improving patient care and quality of life.



**NIPRO MEDICAL CORPORATION**  
200 Crossing Blvd., Bridgewater, NJ 08807  
T: +1 (908) 393-7030  
[www.nipro.com](http://www.nipro.com) | [info@nipromed.com](mailto:info@nipromed.com)

© 2022 Nipro Medical Corp. All Rights Reserved.  
NMC-LBL-0170\_24 March 2022-USA-EN

# OKAY<sup>®</sup> II & GM-30

## Y-Connector & Inflation Device



[www.nipro.com](http://www.nipro.com)



# OKAY® II

Y-Connector

## Control made easy Maximum efficiency with minimal effort

OKAY® II includes both a fixed and hemostatic valve with a large 10 F diameter and one-touch control that facilitates a variety of interventional devices, improving your workflow with minimal effort.




### Single-handedly improving your workflow

- One-touch control for single-hand operation
- “Lock-open” feature to easily clear air and clotting in the system
- Hemostatic valve to minimize blood loss during procedures
- Maximum pressure resistance of 500psi to prevent loss of pressure\*

### Easy to use & ergonomic

- Easy insertion of a variety of multi-device techniques thanks to 10 F inner lumen diameter
- Large grip area for more ergonomic, comfortable hold
- Elongated side port offers a secure grip between fingers

## Valve operation

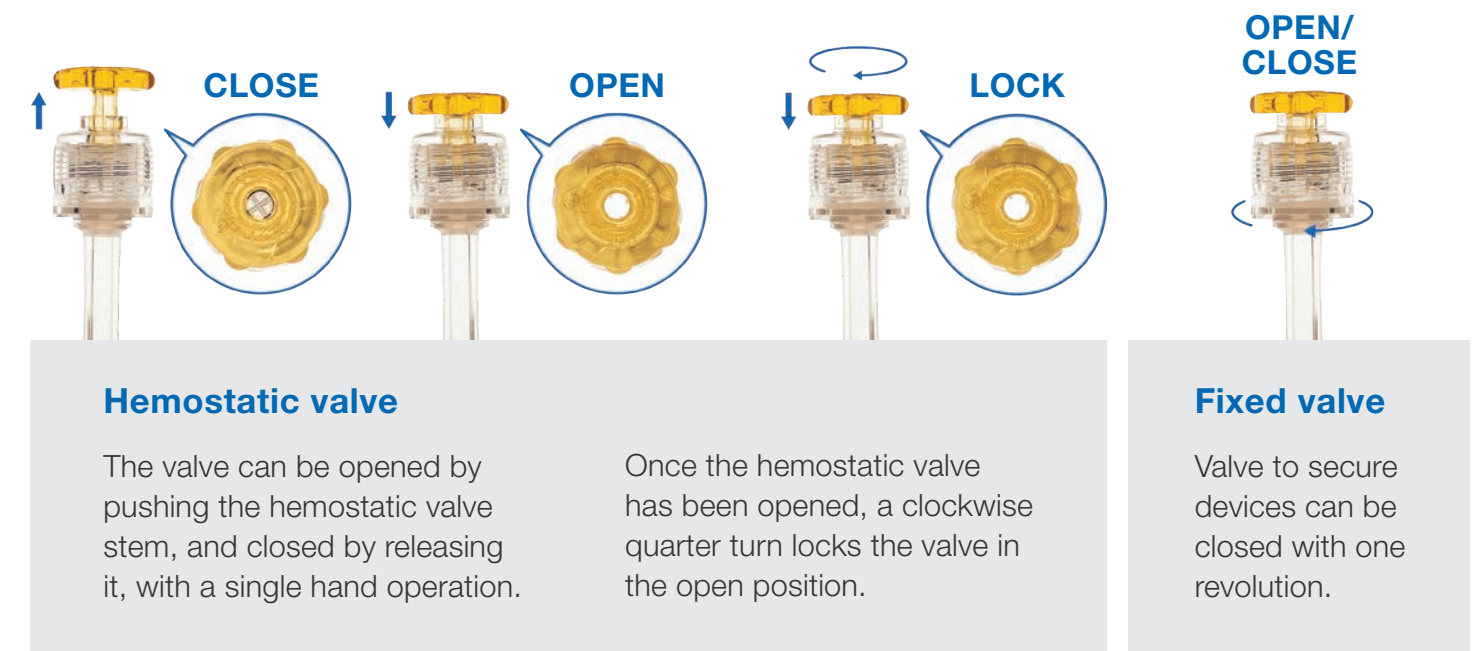

## Accessories

### 1. Extension tube


Length: 22 cm  
Pressure resistance: ≤500psi with fixed valve closed\*

### 2. Inserter

Inner diameter: 0.021"  
Max guidewire size: 0.018"

### 3. Torque device

Compatible guidewire diameters: 0.010"-0.018"

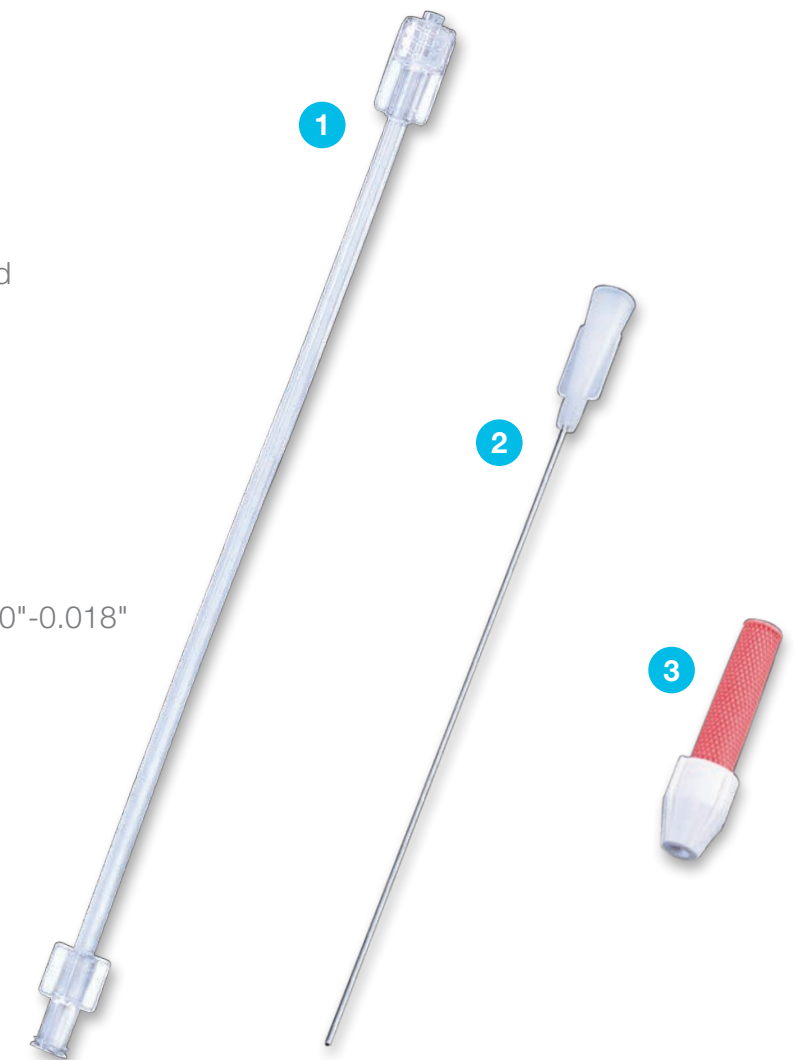

\*Measurements were obtained without any devices inserted. Use of devices (including multiple devices) may result in instances whereby 500psi is not achieved.

# GM-30

Inflation Device

## Control made easy Maximum efficiency with minimal effort

GM-30 has an easy-to-use release lever that requires minimal effort for deflation along with an ergonomically designed handle that allows you to achieve comfortable and confident control with less force and fewer turns.

### Ease of control with minimal effort

- Minimal effort required for inflation and deflation
- Ergonomic handle fits comfortably in the hand, designed for easy rotation while inflating and deflating
- Features a 30 ATM pressure rating with a 25 mL syringe volume

### Achieve maximum performance

- Suitable for many balloon sizes and diameters
- Highly visible red ring on plunger gives clear indication for easy reading of numbers
- Available in two stopcock styles and two tubing lengths
- Kink-resistant tubing



### Ergonomically-designed handle



Less force and fewer turns

Smooth inflation with 1:1 response\*

\*internal testing results on file


### Easy-to-use release lever



Minimal effort required for deflation

Suitable for many balloon sizes and diameters

### Options for various interventional procedures

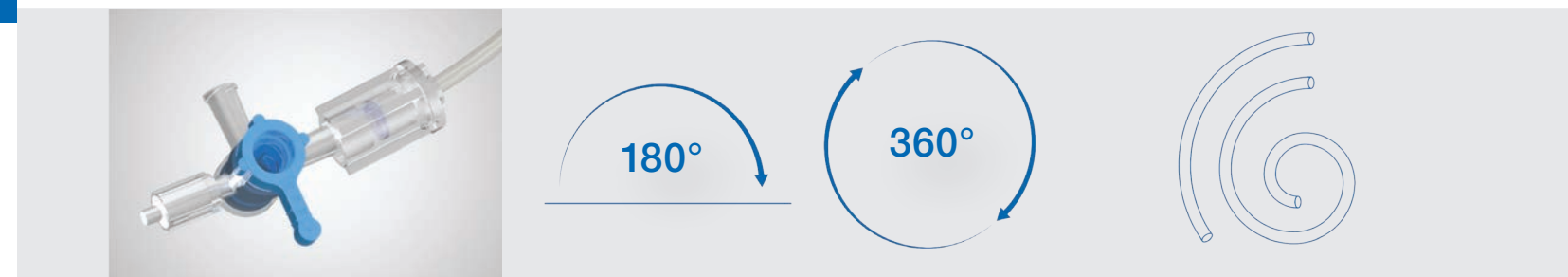

Two stopcock styles:  
3-way and 4-way

Two tubing lengths:  
31 cm and 76 cm

# COMBINATION KIT

OKAY® II & GM-30

Optimize your efficiency in the Cath lab with comfort and ease thanks to these top-performing devices – now available individually or in one convenient kit.

Each combination kit comes with individually-packaged sterile components – offering you the freedom to use both products at the time of procedure or save one sterile component for a later point in time.

## Maximize convenience and minimize waste

- More convenient to bundle
- Helps you be more efficient with inventory
- Reduce waste and save on costs in doing so

## Ordering information

Product code	Product name	Description
YOK0A	OKAY® II Connector Set	Y-Connector, inserter, and torque device
YOK0B	OKAY® II Connector	Y-Connector
YOK0E	OKAY® II Connector Set with extension tube	Y-Connector, inserter, torque device, and extension tube
GM-30N	GM-30 Inflation Device (30atm/25ml)	31 cm tube and 3-way stopcock
GM-30NL	GM-30 Inflation Device (30atm/25ml)	76 cm tube and 3-way stopcock
GM-30NF	GM-30 Inflation Device (30atm/25ml)	31 cm tube and 4-way stopcock
GLGMN-YOK0A	Combination Kit*	<ul style="list-style-type: none"> <li>■ YOK0A (OKAY® II Y-connector, inserter, torque device)</li> <li>■ GM-30N (inflation device, 3-way stopcock, 31 cm)</li> </ul>
GLGMN-YOK0E	Combination Kit*	<ul style="list-style-type: none"> <li>■ YOK0E (OKAY® II Y-connector, inserter, torque device, extension tube)</li> <li>■ GM-30N (inflation device, 3-way stopcock, 31 cm)</li> </ul>
GLGMF-YOK0A	Combination Kit*	<ul style="list-style-type: none"> <li>■ YOK0A (OKAY® II Y-Connector, inserter, torque device)</li> <li>■ GM-30NF (inflation device, 4-way stopcock, 31 cm)</li> </ul>
GLGMF-YOK0E	Combination Kit*	<ul style="list-style-type: none"> <li>■ YOK0E (OKAY® II Y-connector, inserter, torque device, extension tube)</li> <li>■ GM-30NF (inflation device, 4-way stopcock, 31 cm)</li> </ul>

\*Each device comes individually packaged in a sterile pouch