

OKAY[®] II & GM-30

Y-Connector & Inflation Device

YOUR LINK BETWEEN PATIENT AND CARE

Nipro Medical Corporation is a subsidiary of Nipro Corporation Japan, established in 1954 and a world-leader in healthcare with over 30,000 employees around the globe. We are committed to delivering innovative, high quality solutions in the form of renal, medical-surgical, and vascular products. This fundamental objective drives us to develop cutting-edge, value-added products for healthcare professionals that are safe, effective, and easy to use.

As an integral part of a vast global manufacturing network, Nipro Medical Corporation's North America division is responsible for managing the sales, marketing, and business operations of the United States and Canada – unified in our commitment to improving patient care and quality of life.



NIPRO MEDICAL CORPORATION
200 Crossing Blvd., Bridgewater, NJ 08807
T: +1 (908) 393-7030
www.nipro.com | info@nipromed.com

© 2022 Nipro Medical Corp. All Rights Reserved.
NMC-LBL-0170_24 March 2022-USA-EN

OKAY[®] II & GM-30

Y-Connector & Inflation Device



www.nipro.com



OKAY® II

Y-Connector

Control made easy Maximum efficiency with minimal effort

OKAY® II includes both a fixed and hemostatic valve with a large 10 F diameter and one-touch control that facilitates a variety of interventional devices, improving your workflow with minimal effort.




Single-handedly improving your workflow

- One-touch control for single-hand operation
- “Lock-open” feature to easily clear air and clotting in the system
- Hemostatic valve to minimize blood loss during procedures
- Maximum pressure resistance of 500psi to prevent loss of pressure*

Easy to use & ergonomic

- Easy insertion of a variety of multi-device techniques thanks to 10 F inner lumen diameter
- Large grip area for more ergonomic, comfortable hold
- Elongated side port offers a secure grip between fingers

Valve operation

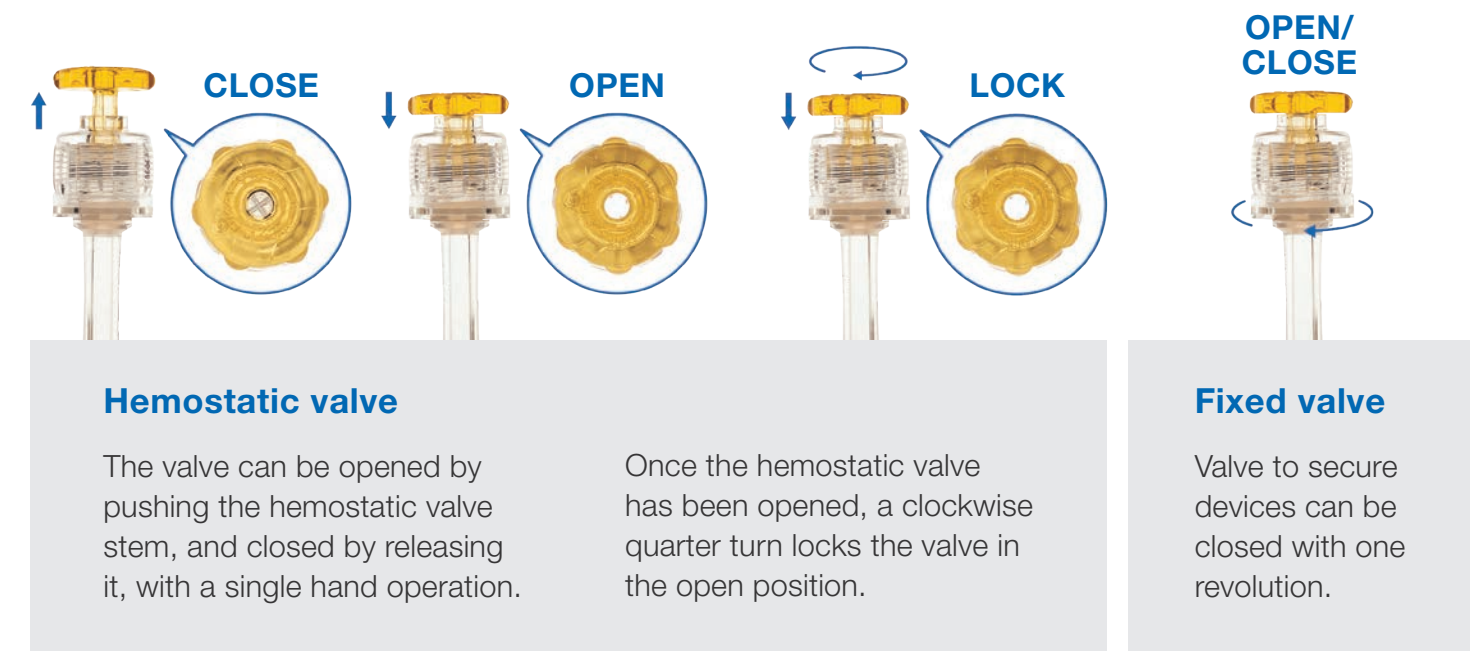

Accessories

1. Extension tube


Length: 22 cm
Pressure resistance: ≤500psi with fixed valve closed*

2. Inserter

Inner diameter: 0.021"
Max guidewire size: 0.018"

3. Torque device

Compatible guidewire diameters: 0.010"-0.018"

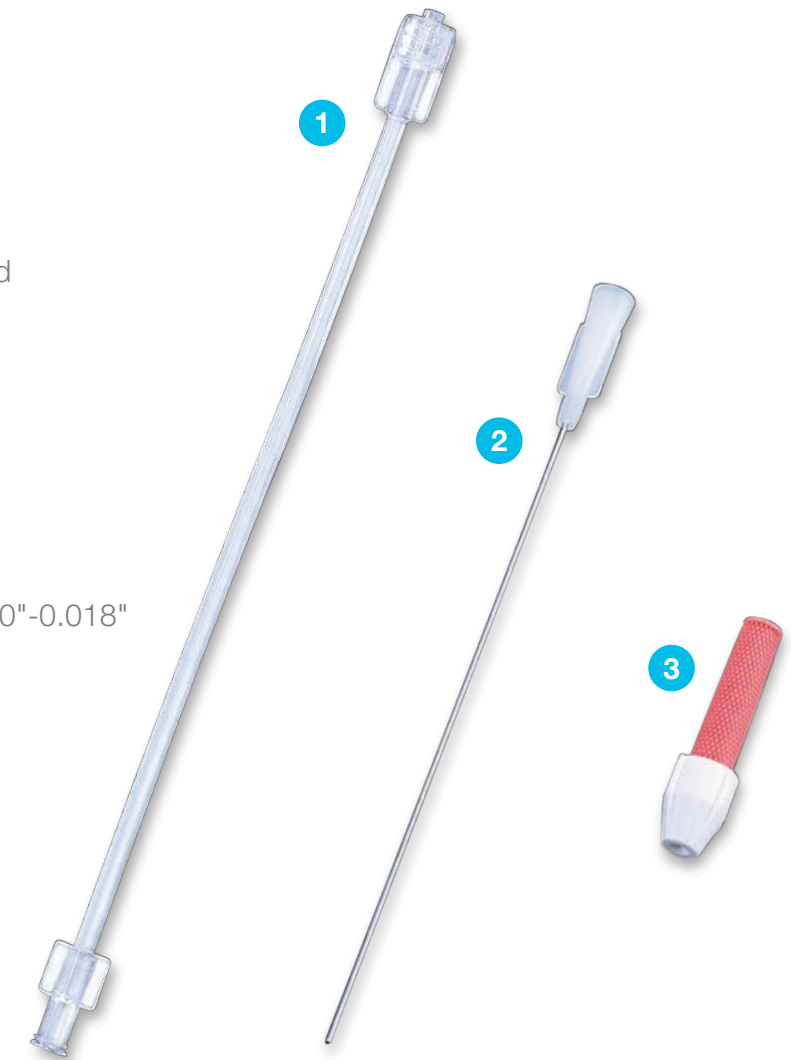

*Measurements were obtained without any devices inserted. Use of devices (including multiple devices) may result in instances whereby 500psi is not achieved.

GM-30

Inflation Device

Control made easy Maximum efficiency with minimal effort

GM-30 has an easy-to-use release lever that requires minimal effort for deflation along with an ergonomically designed handle that allows you to achieve comfortable and confident control with less force and fewer turns.

Ease of control with minimal effort

- Minimal effort required for inflation and deflation
- Ergonomic handle fits comfortably in the hand, designed for easy rotation while inflating and deflating
- Features a 30 ATM pressure rating with a 25 mL syringe volume

Achieve maximum performance

- Suitable for many balloon sizes and diameters
- Highly visible red ring on plunger gives clear indication for easy reading of numbers
- Available in two stopcock styles and two tubing lengths
- Kink-resistant tubing



Ergonomically-designed handle



Less force and fewer turns

Smooth inflation with 1:1 response*

*internal testing results on file


Easy-to-use release lever



Minimal effort required for deflation

Suitable for many balloon sizes and diameters

Options for various interventional procedures

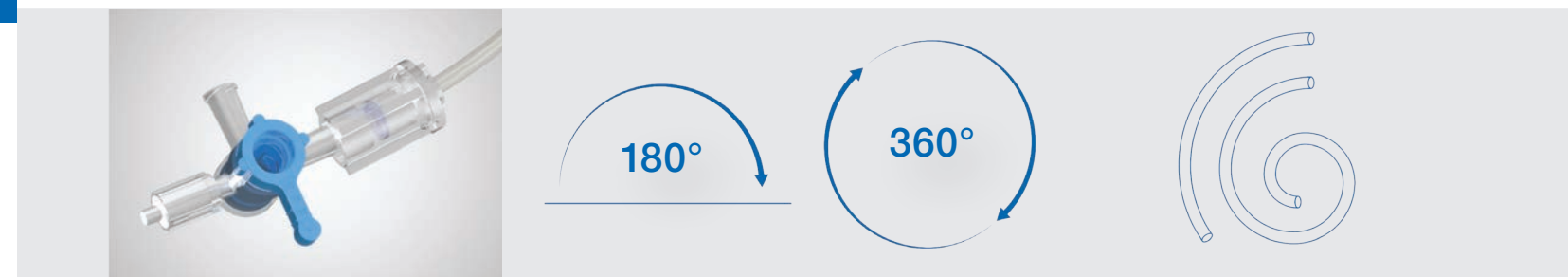

Two stopcock styles:
3-way and 4-way

Two tubing lengths:
31 cm and 76 cm

COMBINATION KIT

OKAY® II & GM-30

Optimize your efficiency in the Cath lab with comfort and ease thanks to these top-performing devices – now available individually or in one convenient kit.

Each combination kit comes with individually-packaged sterile components – offering you the freedom to use both products at the time of procedure or save one sterile component for a later point in time.

Maximize convenience and minimize waste

- More convenient to bundle
- Helps you be more efficient with inventory
- Reduce waste and save on costs in doing so

Ordering information

Product code	Product name	Description
YOK0A	OKAY® II Connector Set	Y-Connector, inserter, and torque device
YOK0B	OKAY® II Connector	Y-Connector
YOK0E	OKAY® II Connector Set with extension tube	Y-Connector, inserter, torque device, and extension tube
GM-30N	GM-30 Inflation Device (30atm/25ml)	31 cm tube and 3-way stopcock
GM-30NL	GM-30 Inflation Device (30atm/25ml)	76 cm tube and 3-way stopcock
GM-30NF	GM-30 Inflation Device (30atm/25ml)	31 cm tube and 4-way stopcock
GLGMN-YOK0A	Combination Kit*	<ul style="list-style-type: none"> ■ YOK0A (OKAY® II Y-connector, inserter, torque device) ■ GM-30N (inflation device, 3-way stopcock, 31 cm)
GLGMN-YOK0E	Combination Kit*	<ul style="list-style-type: none"> ■ YOK0E (OKAY® II Y-connector, inserter, torque device, extension tube) ■ GM-30N (inflation device, 3-way stopcock, 31 cm)
GLGMF-YOK0A	Combination Kit*	<ul style="list-style-type: none"> ■ YOK0A (OKAY® II Y-Connector, inserter, torque device) ■ GM-30NF (inflation device, 4-way stopcock, 31 cm)
GLGMF-YOK0E	Combination Kit*	<ul style="list-style-type: none"> ■ YOK0E (OKAY® II Y-connector, inserter, torque device, extension tube) ■ GM-30NF (inflation device, 4-way stopcock, 31 cm)

*Each device comes individually packaged in a sterile pouch